



KeriSafe

Range of single-use surgical knives & instruments



Guarantees . . .

- consistent guided cutting
- precise control throughout the ligament cutting procedure
- median nerve protection
- palmar arch protection



KeriKnife™
Surgical knife & guide



Mini KeriKnife™
Mini surgical knife

Fully committed to **Hand & Wrist** only


KeriKnife™

Surgical knife & guide


single-use kit, for **carpal tunnel release**

Surgical knife & guide specially developed
for *ultrasound-guided surgery*



Blade design specifically tailored to cut
the transverse carpal ligament, regardless of thickness

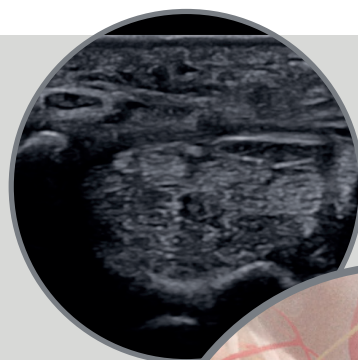

Optimum ultrasound-guided visualisation

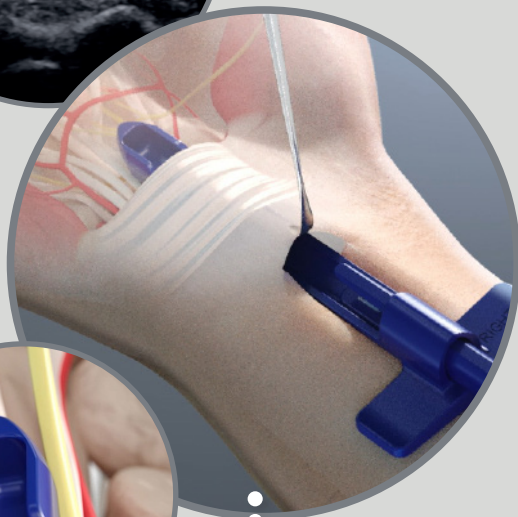

Safety with precisely controlled & consistent
cutting under ultrasound guidance

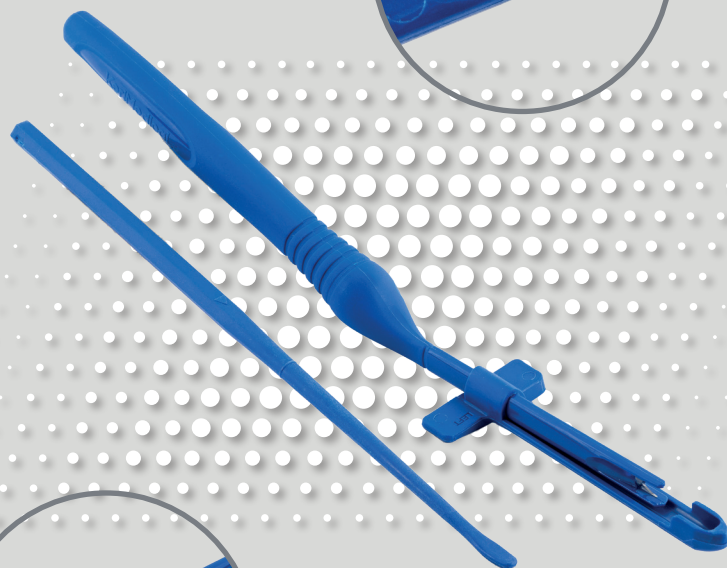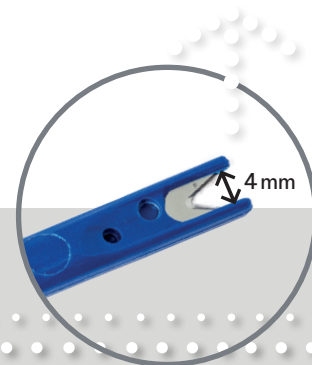

Sterile kit

For preparation and complete
release of the transverse
carpal ligament

Instrument advantages

Inclined cutting guide, in a dedicated right-hand or left-hand format according to the angle



Right-hand guide or left-hand guide

Isolates the median nerve from the cutting path, which ensures it remains protected from the knife

Median nerve protection



Cannula

Guides the knife along its path



Head



Halts the knife, enabling precise control throughout the ligament cutting procedure

Palmar arch protection

Mini Keriknife™

Mini Surgical Knife

single-use kit, for *ultrasound-guided* trigger finger release

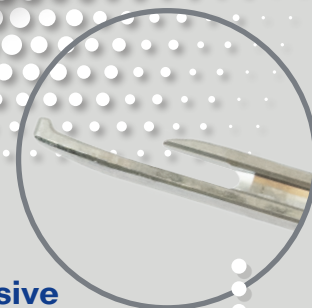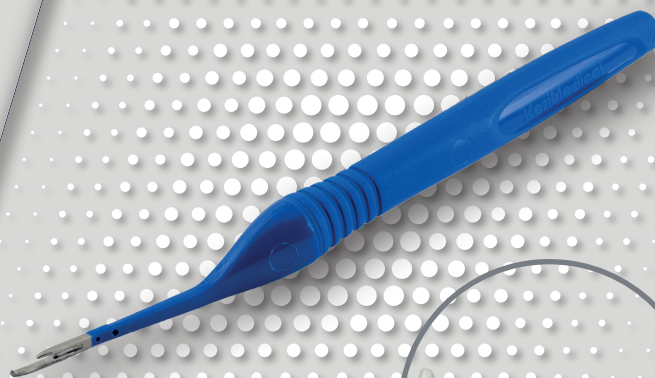

Instrument
ergonomic & minimally invasive

Blade visible with ultrasound guidance

Lower tip: enables only the A1 pulley to be cut and disengages the blade before it reaches the A2 pulley

The advantages OF SINGLE USE



Single-use kit: minimises contamination risk, with full traceability






Sterile and ready to use: the ideal solution for daily operations, with no limits on the number of procedures



An economical solution: reduced sterilisation costs



Convenient sales unit: pack of 5

	Catalogue	Sales unit
Carpal tunnel	Surgical knife, right guide & spatula 524-A90000 	5
	Surgical knife, left guide & spatula 524-A91000 	5
Trigger finger	Mini surgical knife 524-A80000 	5



KERI MEDICAL SA

Chemin du Pré-Fleuri 5
 1228 Plan-les-Ouates - Switzerland
 Tél: +41 58 255 01 30
 Fax: +41 58 255 01 40
 customerservice@kerimedical.com



KERI MEDICAL FRANCE

34, rue Antoine Redier
 74160 Archamps - France
 Tel.: +33(0)4 50 85 36 10
 Fax: +33 (0)4 80 16 05 10
 adv@kerimedical.com

